


DUO DIVISYS®



Bed on fixed feet with MÉDIDOM II* Light Dragon Tree color trim, with long side panels and skirting
*not compatible with side rails



BED FOR COUPLES UP TO 270 KG

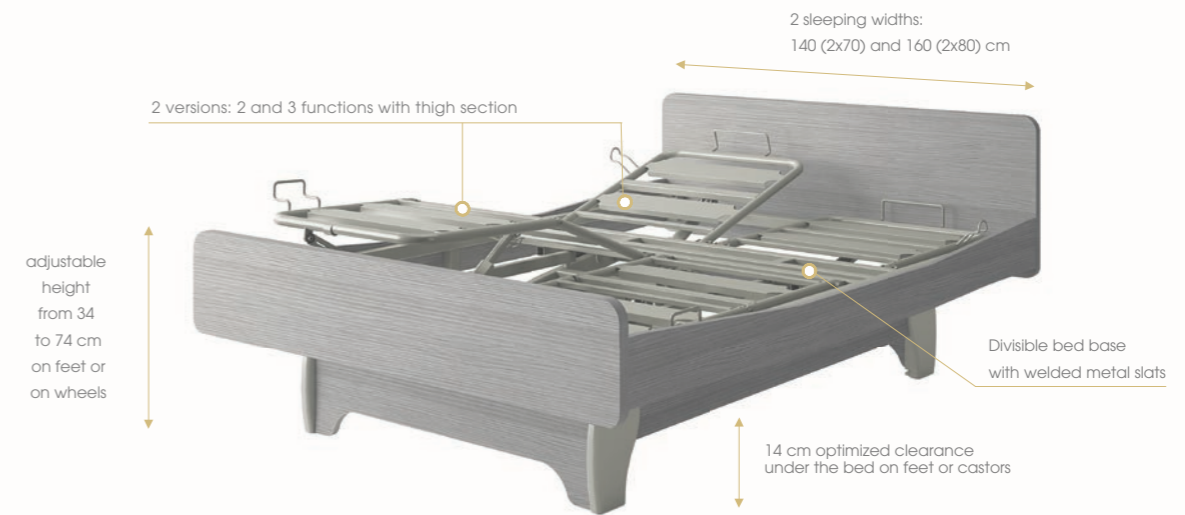

Auzence II board, available in 140 and 160 cm.



Abélia II board, available in 140 and 160 cm.

Bed on fixed feet with MÉDIDOM II* Grey Dragon Tree color trim, with long side panels and skirting.
*incompatible with rails.

3 basic finishes:

On Castors / On castors with skirting / On fixed feet with skirting.



Options:

- ▶ Wooden slatted bed base.
- ▶ Two infrared remote controls.

LIVE TOGETHER FOR AS LONG AS POSSIBLE WHILE RESPECTING EACH OTHER'S INDIVIDUAL NEEDS. REST, READ, EAT OR SLEEP.

PERSONALIZE YOUR HOME FURNITURE. DECK CHAIR, HALF-SEATED OR FLAT POSITION.

DUO DIVISYS BED

- ▶ For couples (combined weight of up to 270 kg), with two separate bed bases for a combined width of 140 cm (2x70) or 160 cm (2x80).

COMFORT AND HYGIENE

- ▶ Divisible bed base with welded metal slats.

ERGONOMICS

- ▶ Adjustable height from 34 to 74 cm.

DESCRIPTION	DIMENSIONS	ADJUSTABLE HEIGHT	BR	LEG REST	LIFTING CAPACITY	CASTORS	STANDARDS
DUO DIVISYS 2 functions electrical	140 x 200 cm	34 cm to 74 cm	70°	12° manual cranks	SWL ⁽¹⁾ 315 kg mattress 30 kg patient 270 kg accessories 15 kg	Ø 100 mm with brake without thread guard	IEC 60601-2-52
	160 x 200 cm						
DUO DIVISYS 3 functions electrical	140 x 200 cm	34 cm to 74 cm	70°	26° thigh sections (thighs) and 16° Feet adjustment by manual cranks	SWL ⁽¹⁾ 315 kg mattress 30 kg patient 270 kg accessories 15 kg	Ø 100 mm with brake without thread guard	IEC 60601-2-52
	160 x 200 cm						

⁽¹⁾ Safe working load

- 2 accessory holders (lengthwise)
- 2 IV stand holders (at the head of the bed)
- 2 lifting pole (at the head of the bed)

Handset, electric and mechanical parts: 5-year warranty*

Boards and accessories: 1-year warranty*

* Under normal conditions of use and excluding wearing parts - please see instructions of use

Photos not contractually binding

FP000007/01.2018

WINNOCARE Group - 4, Le Pas du Château - 85670 Saint-Paul-Mont-Penit - France

Tél : +33 (0)2 51 98 55 64 - Fax : +33 (0)2 51 98 59 07 - E-mail : info@medicatlantique.fr - www.winnocare.fr

In the interest of constantly improving its products, the WINNOCARE Group reserves the rights to modify the technical specifications of any product without prior notice

Class I medical device in accordance with European Directive 93/42/EEC

